

Meridium. Reclaim your way.

The only true microprocessor-controlled foot



What is new?



Meridium is now also compatible with iOS
Cockpit App available for Android and iOS devices



Improved battery life
Up to 44 % longer runtime



Quality improvements
New overheating mode

Technical data


Suitable for	MG 2-3
Max. body weight	125 kg / 125 lbs
Sizes	24-29 cm
Color options footshell	translucent, light brown, beige
Approved knee joints	3C98* C-Leg, 3B5* Genium X3, 3B1* Genium

More information about the update at a glance


- Article number stays unchanged (1B1-2=*)
- No new certification needed
- Update of the M-Soft to version 1.4 is necessary (4X154=V1.4)
- Remote control (4X350) is no longer needed for iOS compatible Meridiums

What makes Meridium unique?


Unique adaptation to the human gait cycle

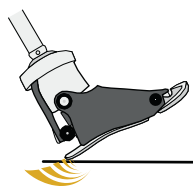

▶ Quick and stable ground contact

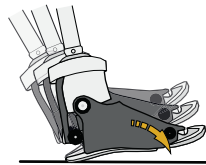



▶ Adaptation of the ankle allows a more natural rollover

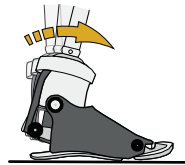

▶ Longer ground contact thanks to flexible toe plate



▶ Reduced risk of tripping due to increased ground clearance in the swing phase

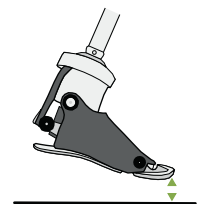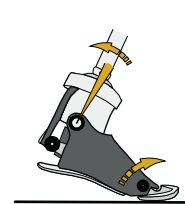

Natural walking thanks to real-time adaptation and movability up to the toes



Automatic heel height adjustment up to 5 cm



Safer down stairs walking with full step contact



Stable standing even on slopes due to an innovative stance function